

IPL SHR Multi Treatment System



JDS-100

— *Machine Introduction* —



- Hair Removal
- Skin Rejuvenation
- Acne Removal



— Working Principle —

IPL is a non-invasive skin treatment that emits powerful bursts of light energy to treat sun damaged skin, wrinkles, and unwanted hair growth, etc. without damaging the surface of your skin.

IPL systems use hand held flash lamps to emit micro pulses of light across a broad band of visible wavelengths. By using gentle pulses of light momentarily to create heat, IPL can permanently disable hair follicles. and prevent future growth. Similarly, IPL also achieves excellent long term results when treating thread veins, skin pigmentation, acne and other skin complaints.

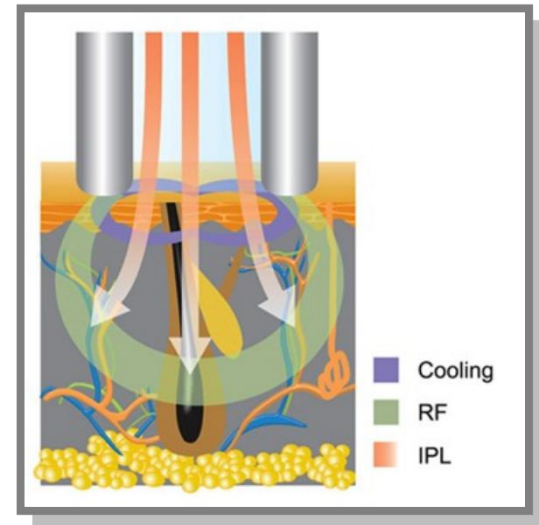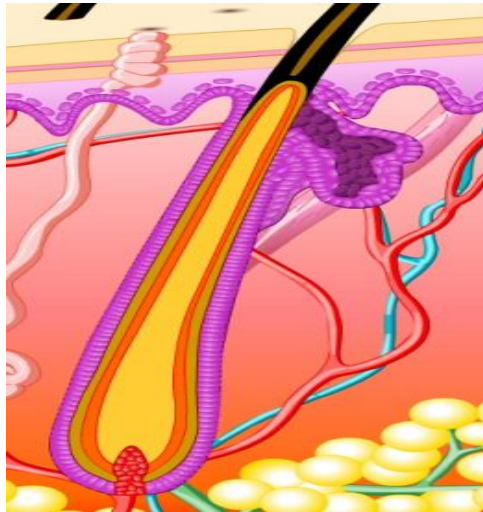

—Handle Introduction—

Big Spot 15*50mm



Switch



➤ **Standard:** 430/530/640-1200nm

➤ **Optional:** 590/690/750-1200nm

— Before and After —



Female | Age 26
Before

After
5 treatments

Male | Age 39
Before

After
4 treatments



— Machine Parameters —

Power Supply	2000W
Screen	10.4 Inch TFT Screen
Wavelength	IPL(430/530/640-1200nm),SHR(650-950nm)
Energy	1-26J/cm ²
Spot Size	15*50mm
Pluse Amount	1-6
Pluse	1-9.9ms
Cooling System	Air Cooling + Water Cooling + TEC
Machine Size	52*38*37
Net Weight	35Kgs

—About JDS—

JDS Laser has been servicing distributors and customers across different countries for 10 years. A combination of expertise, passion and a sense for great customer care has led the company to grow from 1 location to 2 while expanding from Beijing to a second one located in Handan, as well as a soon-to-be 3rd location (Lisbon, Portugal) and 4th location (Toronto, Canada). JDS Laser rapidly became one of the best sell Beauty Equipment Supplier in the region, earning itself multiple awards of the Year.

—Our Mission—

Our company was founded on the concept of selling Beauty Equipment that perfectly blends our selling skills in the Beauty Equipment industry with next-level customer care. We understand the industry and know how it continuously changes, which is why we work hard to stay on top of our competitors by making sure our staff is up to date with the latest information, services and training. We always provide honest recommendations and a multitude of dependable products.



JDS Beauty[®]
Aesthetic ● Technology ● Innovation

Focus on
Quality and Service