

# ICE Diode Laser Hair Removal System



**JDS-808**

# — Machine Introduction —

- Hair Removal
- Skin Rejuvenation




6 Skin Types



## — Working Principle —

Diode laser hair removal system, adopts the light selective absorption theory, heating the hair shaft and hair follicle to destroy the hair follicle and oxygen organization around hair follicle, thereby to get the permanent hair removal purpose.

The patent pulse width to energy 1:1 working mode, makes energy penetrating more effective, less thermal destroy to surrounding tissues, meanwhile, the stable anti-ice cooling technology, protects skin from burny finally getting a wonderful result.

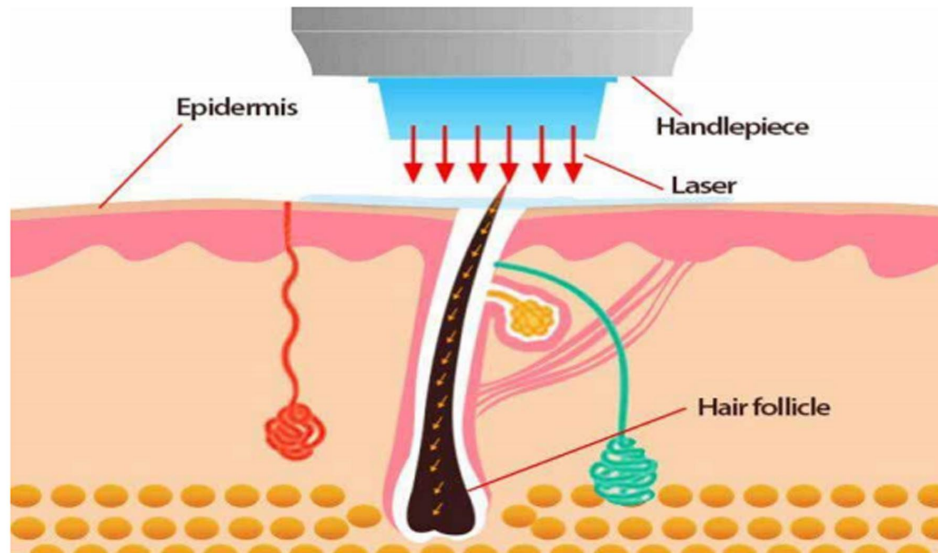

# — Handle Introduction —

▶ *Long lifespan*  
USA Imported Coherent Laser Bar



▶ *More convenient*  
Handle with Display Screen

## — Before and After —



Female | Age 26  
Before

---

After  
5 treatments

Male | Age 39  
Before

---

After  
4 treatments



## — Machine Parameters —

<b>Laser Type</b>	<b>808nm/808+755+1064nm</b>
<b>Interface</b>	<b>10.4 inch color LCD touch screen</b>
<b>Spot Size</b>	<b>13*16mm/13*22mm/13*30mm</b>
<b>Handle Power</b>	<b>1000W/1200W/1600W</b>
<b>Frequency</b>	<b>1-10Hz adjustable</b>
<b>Pulse Width</b>	<b>10-300ms</b>
<b>Cooling System</b>	<b>TEC+Water+Wind cooling</b>
<b>Voltage</b>	<b>AC220V/110V 50-60Hz</b>
<b>Net Weight</b>	<b>30kg</b>





JDS 捷创达胜科技

**JDS Beauty**<sup>®</sup>  
Aesthetic ● Technology ● Innovation

Focus on  
**Quality and Service**